



PERIPHERAL INTERVENTION

DRAINAGE



Thoracentesis



Biliary Drainage



Nephrostomy



Paracentesis



Abscess Drainage



PERIPHERAL INTERVENTION

COMPLETE DRAINAGE SOLUTIONS

Merit Medical offers an integrated suite of products designed to support your drainage procedures every step of the way. We listen to you to **UNDERSTAND** your needs and **INNOVATE** superior products that overcome your toughest clinical challenges. We consistently provide the quality products you need, when you need them, to help you **DELIVER** efficient, patient-focused, and cost-effective care.



MAK-NV™
Introducer System

ACCESS

Ideal for:



Abscess Drainage




Biliary Drainage



Nephrostomy



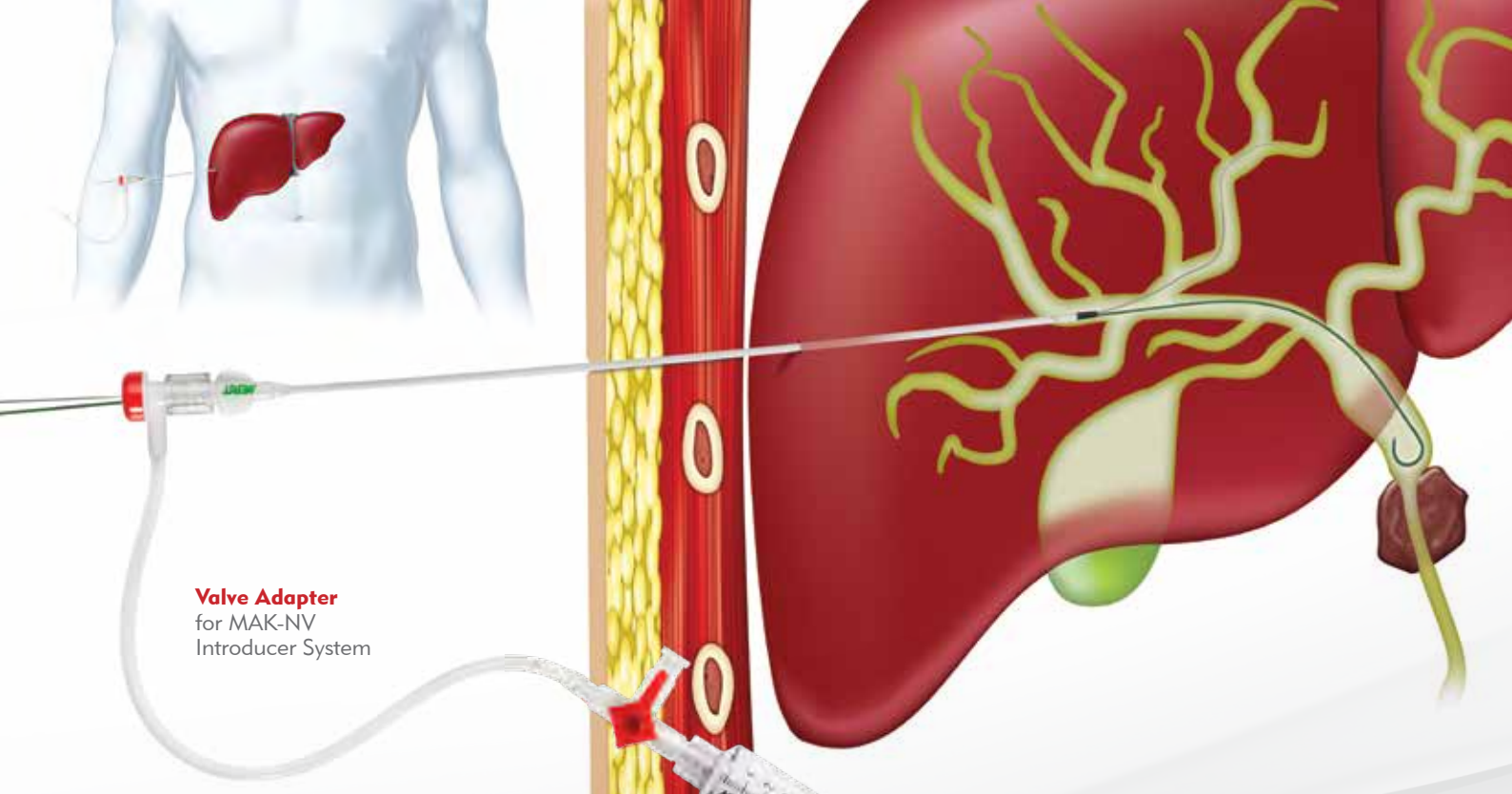

Valve Adapter
for MAK-NV
Introducer System

START WITH QUALITY

Our MAK-NV (Mini Access Kit, Non-Vascular) Introducer System provides quick access and easy **VISUALISATION** into the drainage area with an echo-enhanced needle and radiopaque marker tip on the introducer. Superior tip design **DECREASES BUCKLING** and introducer peel-back for smooth insertion and **EXCELLENT TRANSITIONS**. The Valve Adapter is designed to be back-loaded over a guide wire and attached to the MAK-NV Introducer System to facilitate aspiration or infusion.



Introducer Needle
with Chiba Stylet



Introducer Needle
with Trocar Stylet

DRAINAGE CATHETERS

ReSolve® Biliary
Locking Drainage
Catheter

Two drainage hole patterns to
accommodate multiple lesion locations
and varying patient anatomy

2 cm

5 cm

Unique locking
mechanism to
prevent movement


Ideal for:





Biliary Drainage

Radiopaque marker
band helps with fast and
accurate positioning

Repositioning tool to
unlock and
reposition catheter

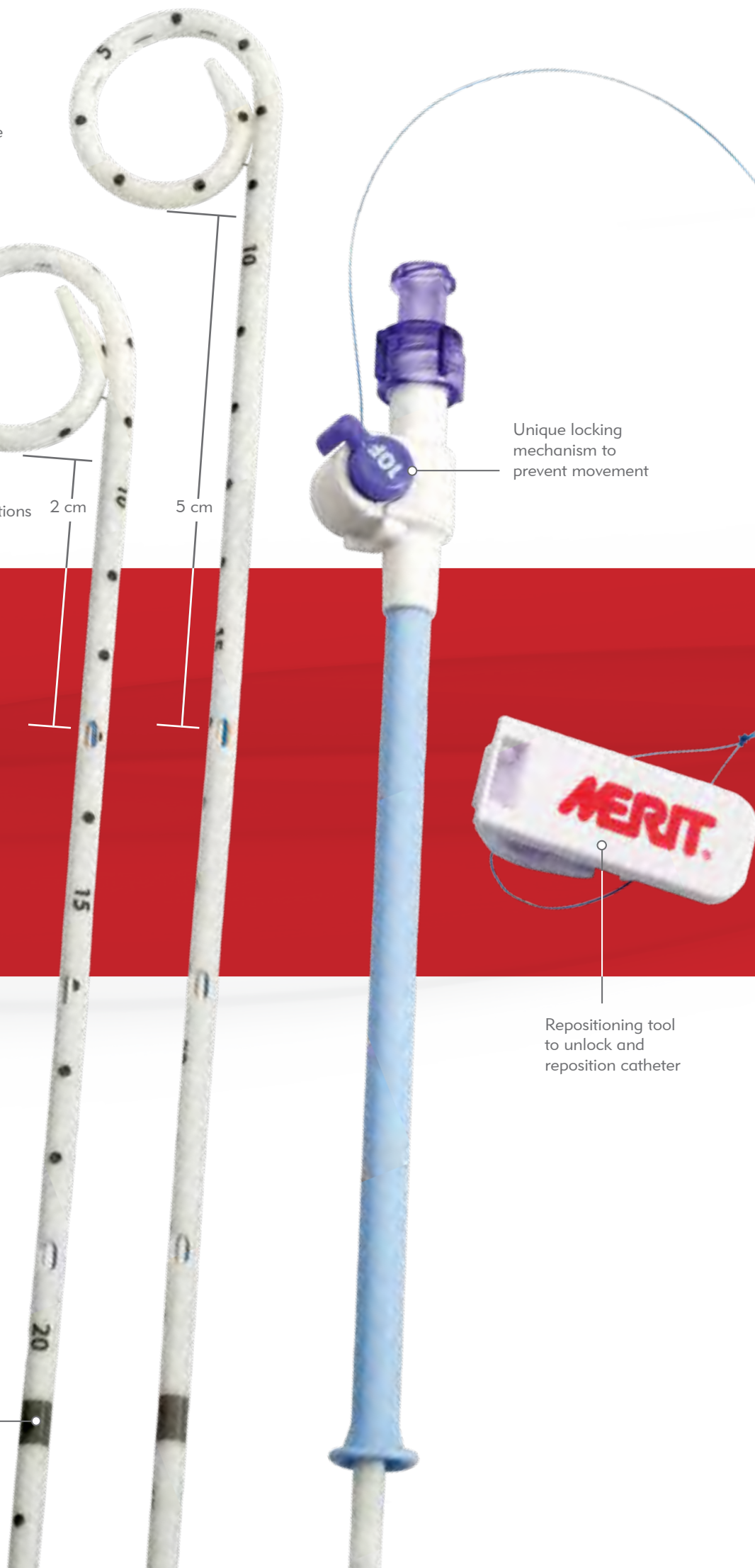


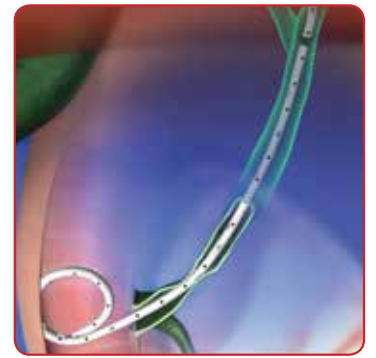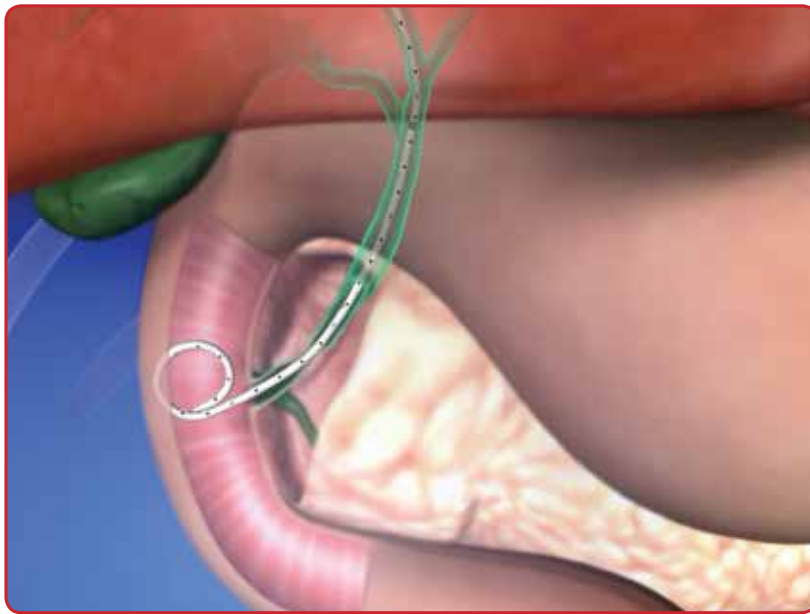

The different hole patterns of the ReSolve Biliary Locking Drainage Catheter provide targeted drainage of both distal lesions (top right illustration, RBDC option) and proximal lesions (bottom right illustration, RBC option).

POWERFUL DRAINAGE

Our ReSolve Biliary Locking Drainage Catheter is designed for **FAST, UNIMPEDED DRAINAGE** and offers up to 4X the drainage area, with multiple **LARGE HOLES** and different size configurations. Additional benefits include depth markers and a radiopaque marker band for **FAST AND ACCURATE POSITIONING**, an atraumatic pigtail designed to reduce the potential for damaging the duodenum, and its polyurethane material allows for long-term indwelling time.

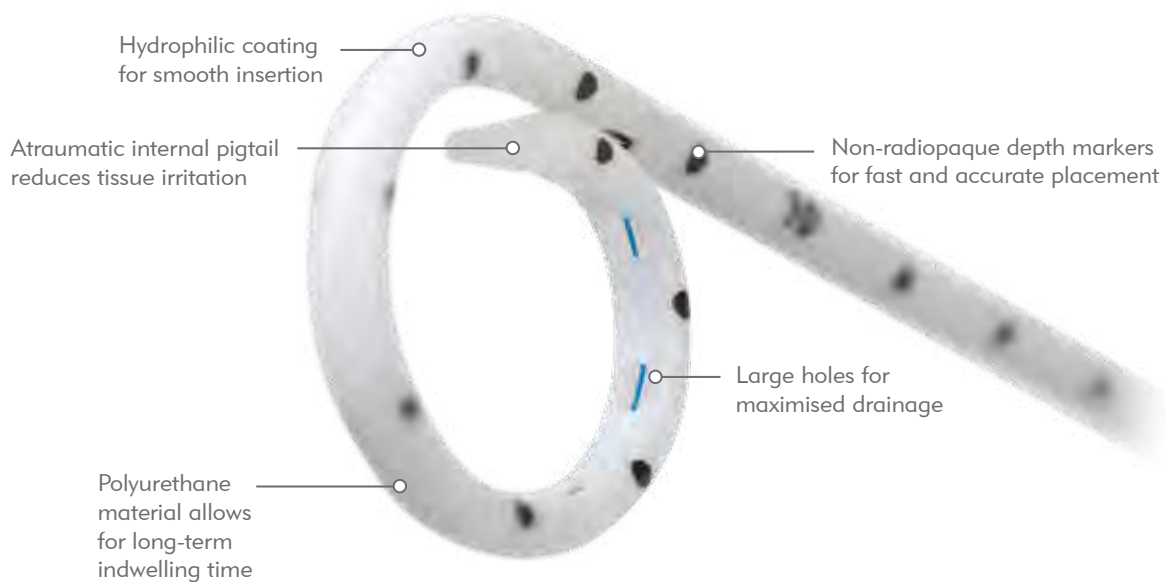

ReSolve® Locking
Drainage Catheter

Unique locking mechanism
to prevent movement

Repositioning tool
to unlock and
reposition catheter

DRAINAGE CATHETERS

Non-radiopaque depth markers for
fast and accurate placement

Ideal for:



Abscess Drainage



Nephrostomy

Atraumatic internal
pigtail reduces
tissue irritation

Tapered tip and
hydrophilic coating
for smooth introduction
and reduced
patient discomfort

Radiopaque marker band helps with
fast and accurate positioning



ReSolve® Non-Locking
Drainage Catheter with J Tip

ReSolve® Non-Locking
Drainage Catheter with Straight Tip

Our ReSolve locking and non-locking drainage catheters are designed to provide **MAXIMUM** drainage, lessen patient discomfort, and promote **MORE ACCURATE PLACEMENT**. The tip of the locking drainage catheter is located on the inside of the pigtail curve to **PREVENT** irritation of the surrounding tissue. The non-locking drainage catheter comes in a J or straight tip.

CENTESIS CATHETERS

Ideal for:



Paracentesis



Thoracentesis

Valved One-Step™
Centesis Catheter
with Slip Fit





TREAT EFFECTIVELY

Our One-Step Centesis Catheter is designed for **MAXIMUM FLOW** with four large drainage holes that spiral around the catheter. The low-profile, tapered tip and **SMOOTH** transitions decrease resistance and patient discomfort. The catheter's echo-enhanced needle increases **VISIBILITY** for more accurate placement. The Valved One-Step features a **SELF-CLOSING** valve that reduces the risk of pneumothorax and fluid leakage.

One-Step™
Centesis Catheter
with Fixed Luer



One-Step™
Centesis Catheter
with Slip Fit



DRAINAGE
ACCESSORIES



Marquis®
Large Bore Stopcock



VacLok®
Vacuum Pressure
Syringe



48" Large Bore Tubing
with pinch clamp and
non-vented spike

Merit Drainage Depot™
Drainage Bag



MORE THAN CATHETERS

Our extensive drainage portfolio also includes ancillary products that make drainage more **EFFECTIVE**, safe and timesaving. These products include **NEEDLELESS** catheter fixation devices, vacuum pressure syringes, large bore stopcocks with **LARGE LUMENS** designed to increase flow rates, drainage tubing, and drainage bags with internal anti-reflux valves and adjustable-length tubing.



Revolution™
Catheter Securement
Device



StayFIX®
Fixation
Device

CUSTOM SOLUTIONS

Ideal for:



Abscess Drainage



Nephrostomy

ReSolve®

Locking Drainage
Catheter Tray



OPTIMISE YOUR DRAINAGE PROCEDURE

Trays, kits and packs are **EFFICIENT**, timesaving and cost-effective. Merit Medical offers numerous pre-arranged options, or we will work with you to create an assortment that **UNIQUELY** meets your needs and provides **EVERYTHING** you require to perform the drainage procedures you do most. Contact a Merit Medical Representative to learn more.



Understand. Innovate. Deliver.™



European Corporate Headquarters, Maastricht, The Netherlands

MeritEMEA.com

Merit Medical Europe, Middle East & Africa

Amerikalaan 42, 6199 AE Maastricht-Airport, Netherlands
+31 43 358 82 22

Merit Medical Ireland Ltd.

Parkmore Business Park West, Galway, Ireland
+35 39 170 37 33

Merit Medical France

136 Ave. Joseph Kessel, 78960 Voisins Le Bretonneux, France
+33 (0)1 30 57 24 24

Merit Medical Systems, Inc.

1600 West Merit Parkway, South Jordan, UT, 84095 USA
+1 801 253 1600

Asia, Pacific Rim, Latin America, Caribbean & Canada

+1 801 208 4438

Austria

0800 295 374

Belgium

0800 72 906 (Dutch)
0800 73 172 (Français)

Denmark

80 88 00 24

Finland

0800 770 586

France

0800 91 60 30

Germany

0800 182 0871

Ireland (Republic)

1800 553 163

Italy

800 897 005

Luxembourg

8002 25 22

Netherlands

0800 022 81 84

Norway

800 11629

Portugal

308 801 034

Russia

+7 495 221 89 02

Sweden

020 792 445

Switzerland

+41 225180252
(Deutsch)

+41 225948000
(Français)

+41 225180035
(Italiano)

UK

0800 973 115