

MERIT
SURECROSS[®]
SUPPORT CATHETER



MERITMEDICAL[®]

| trackability, pushability, visibility

➤ trackability, pushability, visibility

Crossing tight, difficult lesions is often challenging and time consuming. Having a support catheter that offers trackability, pushability and visibility gives you the confidence needed to cross occlusions, even in your most challenging case.

Product Features

Translucent Shaft for visualization of blood to verify placement in vessel lumen

3 Platinum Marker Bands to confirm catheter position and lesion sizing

- First marker band 2 mm from tip
- 0.014" and 0.018" second and third marker bands are 15 mm apart
- 0.035" second and third marker bands are 50 mm apart

40 cm Distal Hydrophilic Coating provides smooth tracking through tortuous anatomy


Low Profile Tapered Tip facilitates crossing difficult lesions

Maximum Infusion Pressure 300psi

Guide Wire Compatibility 0.014", 0.018", 0.035"

Multiple Lengths 65 cm, 90 cm, 135 cm, 150 cm



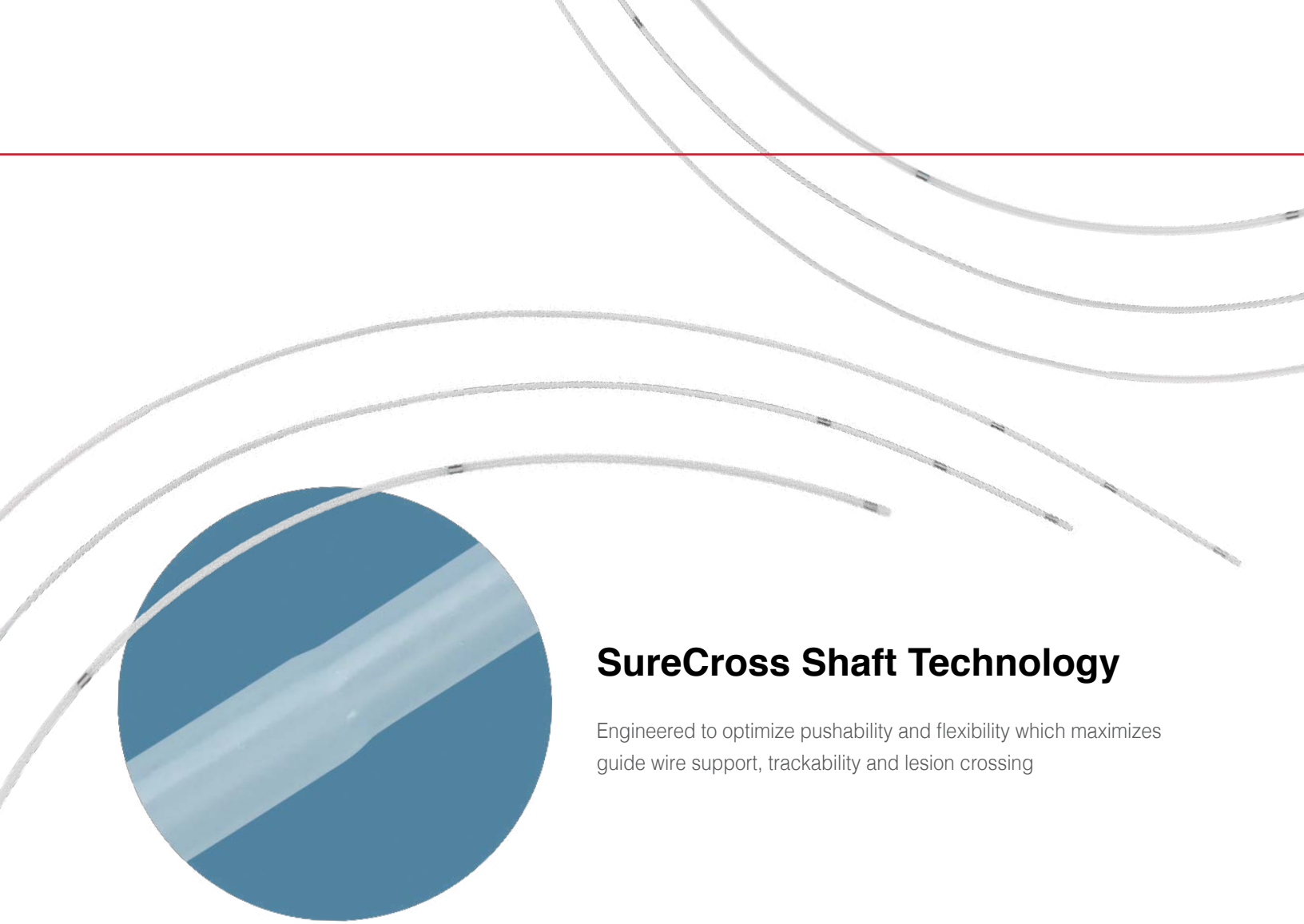

SureCross Shaft Technology

Engineered to optimize pushability and flexibility which maximizes guide wire support, trackability and lesion crossing

SureCross Tip Design

Short tip distance to marker band (2 mm) maximizes accuracy of tip location and enhances tip radiopacity

Internal platinum distal marker band provides a seamless transition zone

Marker band improves critical tip visualization

Robust tapered tip improves crossing and durability


Excellent transition from catheter tip to guide wire

SureCross Tip Design



Tip to Guide Wire Comparison

The low profile, seamless transition from catheter to guide wire provides exceptional crossing support

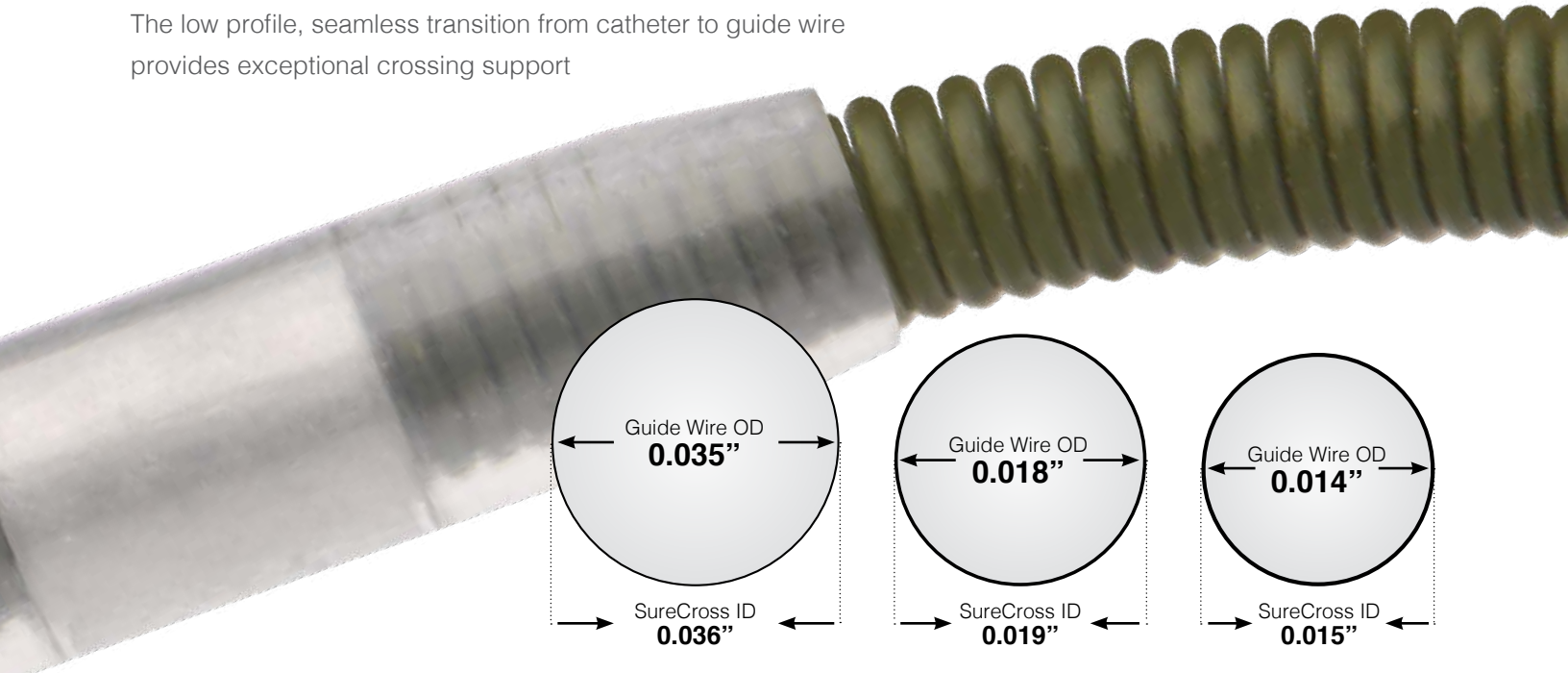

Merit SureCross SUPPORT CATHETER

BowTie™ Guide Wire Insertion Device included in 0.014" and 0.018" configurations



SureCross Catalog #	Guide Wire Compatibility	Minimum Guiding Catheter*	Minimum Sheath Introducer*	Working Length	Number of Marker Bands	2nd and 3rd Marker Band Spacing	Hydrophilic Coating Length	Strain Relief Color	Crossing Profile
SC1465	0.014"	5F	4F	65 cm	3	15 mm	40 cm	Blue	1.5F (3.8 cm)
SC1490	0.014"	5F	4F	90 cm	3	15 mm	40 cm	Blue	1.5F (3.8 cm)
SC14135	0.014"	5F	4F	135 cm	3	15 mm	40 cm	Blue	1.5F (3.8 cm)
SC14150	0.014"	5F	4F	150 cm	3	15 mm	40 cm	Blue	1.5F (3.8 cm)
SC1865	0.018"	5F	4F	65 cm	3	15 mm	40 cm	Purple	1.7F (4.3 cm)
SC1890	0.018"	5F	4F	90 cm	3	15 mm	40 cm	Purple	1.7F (4.3 cm)
SC18135	0.018"	5F	4F	135 cm	3	15 mm	40 cm	Purple	1.7F (4.3 cm)
SC18150	0.018"	5F	4F	150 cm	3	15 mm	40 cm	Purple	1.7F (4.3 cm)
SC3565	0.035"	6F	5F	65 cm	3	50 mm	40 cm	Green	3.1F (7.8 cm)
SC3590	0.035"	6F	5F	90 cm	3	50 mm	40 cm	Green	3.1F (7.8 cm)
SC35135	0.035"	6F	5F	135 cm	3	50 mm	40 cm	Green	3.1F (7.8 cm)
SC35150	0.035"	6F	5F	150 cm	3	50 mm	40 cm	Green	3.1F (7.8 cm)

Packaged 5 per box.

Merit Medical Systems, Inc.
1600 West Merit Parkway
South Jordan, Utah 84095
1-801-253-1600
1-800-35-MERIT
www.merit.com

Merit Medical Europe, Middle East,
and Africa (EMEA)
Amerikalaan 42, 6199 AE Maastricht-Airport
The Netherlands
Tel: +31 43 358 82 22

Merit Medical Ireland Ltd.
Parkmore Business Park West
Galway, Ireland
+353 (0) 91 703 733

Distribuido por:

ventas@cardioparts.com • www.cardioparts.com
Tels. (55) 9183-3933 (81) 8357-0779 (81) 8357-1903